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# STANDING FRAME (FROM 5 TO 12 YEARS) - T.P.1250

It allows to stand safely, stable and independent.

The act of standing helps to stretch the muscles and tendons to prevent contractures, improves the circulation, stretches the abdominal cavity so the respiratory and digestive processes work easily and help in the evacuation of the rectal tract.

Tray adjustable in height and proximity.

Set of hip stabilisers and pelvis waist buffer, adjustable in both width and proximity.

Adjustable leg buffer.

Non-slip feet rest with adjustable feet guards and straps.

Metal frame, with chrome and epoxy finish. All components are padded fireproof upholstery.

## OPTIONAL:

ADJUSTABLE PELVIC SUPPORT – T.P.1250/A

5


USER MANUAL

MANUFACTURER: Ortotecsa,s.l. C/Peña de Francia,8 – 28500 Arganda del Rey (Madrid), España.  
Tif.: +34918719132 / Fax: +34918719323 – email: [info@ortotecsa.com](mailto:info@ortotecsa.com), [www.ortotecsa.com](http://www.ortotecsa.com)

Last revision 10/10/2014

## 5. ASSEMBLY:

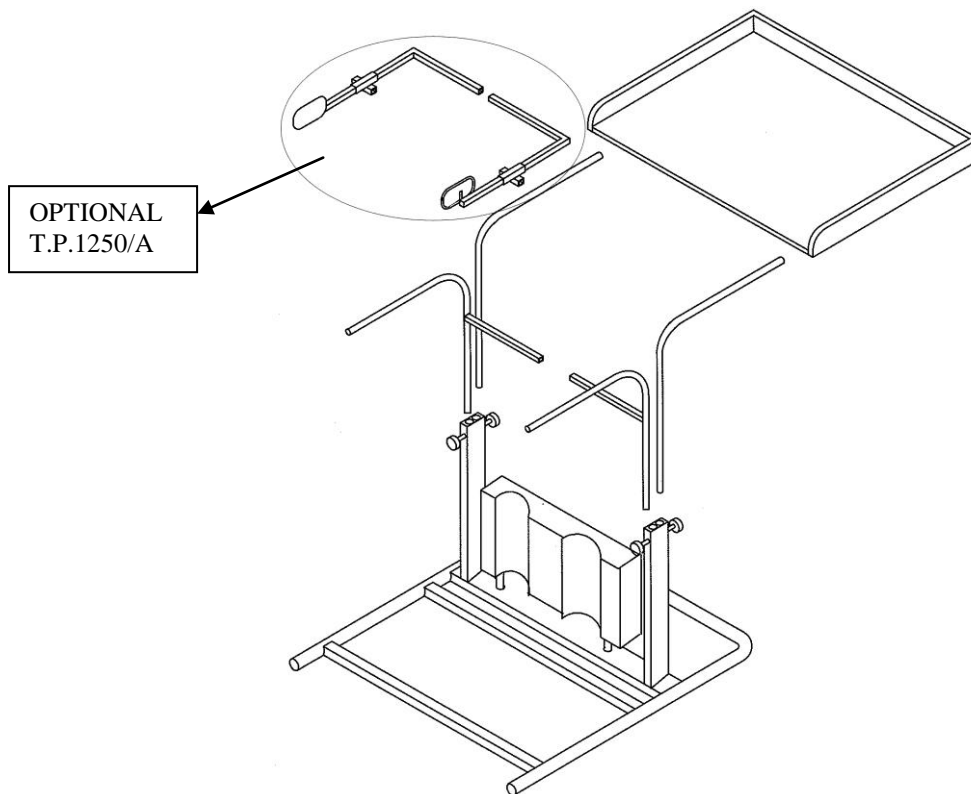Measures: 70x70x50-110 cms.

## 6. HOW TO USE:



Previously graduate all the standing accessories to the right size of the user (SEE point 3: setting accessories). Once the adjustments have been made, you can now use the standing.

Access to the standing by folding down the top of the pelvic waist buffer in the direction shown according to the arrow (SEE Detail \*), place your feet in the feet guards and close the top of the pelvic waist with the fist.

Detail \*

7. ACCESORIES:

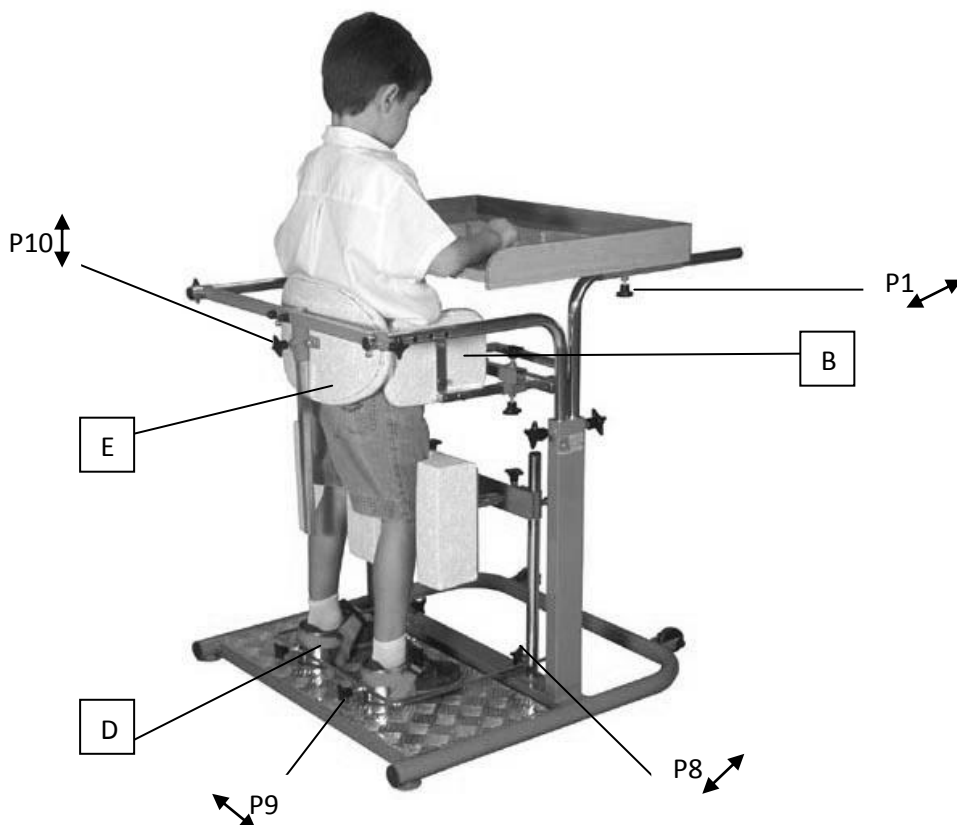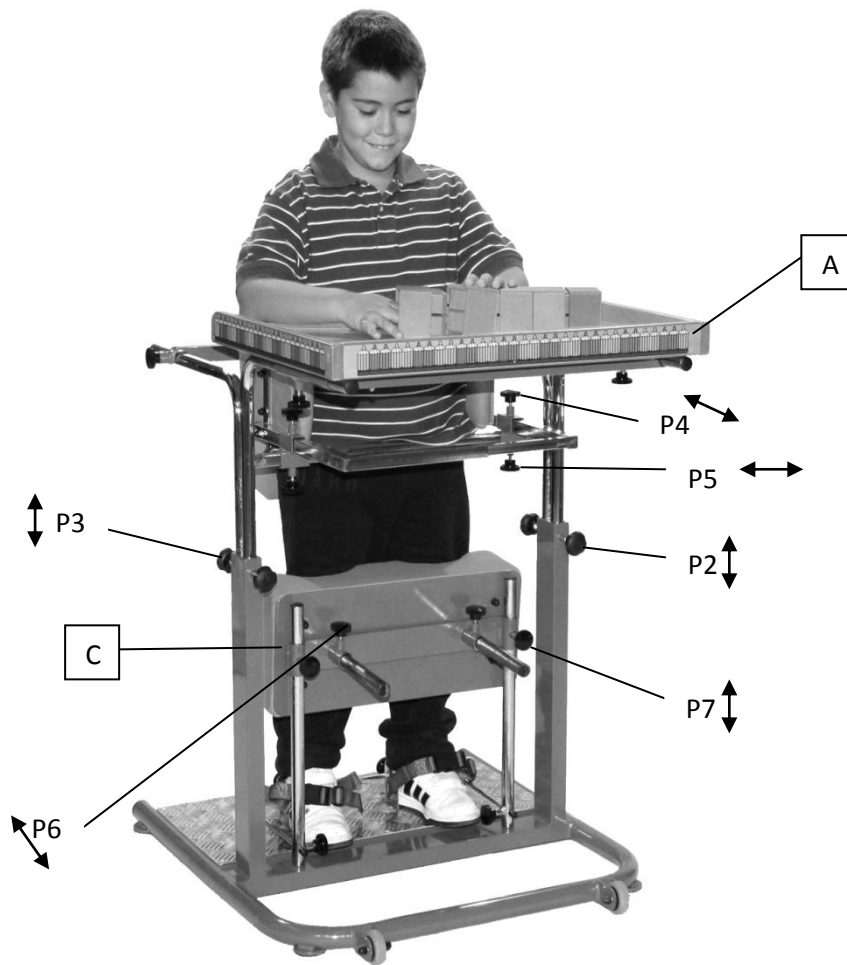- **Graduation tray: (A)**  
**Height:** loosen the knobs P2, slide in the direction of the arrow F2 to the desired position. Set the tray by tightening the knobs P2 again.  
**Depth:** loosen the knobs P1, slide in the direction of the arrow F1 to the desired position. Set the tray by tightening the knobs P1 again.
  
- **Graduation hip controllers or adjustable pelvic support: (B) Optional!**  
**Height:** loosen the knobs P3, slide in the direction of the arrow F2 to the desired position. Set the controller by tightening the knobs P3 again.  
**Depth:** loosen the knobs P4, slide in the direction of the arrow F3 to the desired position. Set the controller by tightening the knobs P4 again.  
**Width:** loosen the knobs P5, slide in the direction of the arrow F4 to the desired position. Set the controller by tightening the knobs P5 again.
  
- **Graduation top legs: (C)**  
**Height:** loosen the knobs P7, slide in the direction of the arrow F6 to the desired position. Set the top by tightening the knobs P7 again.  
**Depth:** loosen the knobs P6, slide in the direction of the arrow F5 to the desired position. Set the top by tightening the knobs P6 again.
  
- **Graduation feet guards feet: (D)**  
**Depth:** loosen the knobs P8, slide in the direction of the arrow F7 to the desired position. Set the feet guards by tightening the knobs P8 again.  
**Width:** loosen the knobs P9, slide in the direction of the arrow F8 to the desired position. Set the feet guards by tightening the knobs P9 again.
  
- **Graduation top pelvic waist (backrest): (E)**  
**Width:** loosen the knob P10, slide in the direction of the arrow to the desired position. Set the backrest by tightening the knob P10 again.

## 8. WARRANTY:

ORTOTECSA, S.L. will provide a warranty for two years from the date of purchase against manufacturing defects. Our warranty includes repair, replacement, or exchange of the product and / or components including labor in our installations.

To have effect this guarantee, is required the purchase invoice.

To have effect this guarantee, is required the purchase invoice.

### **THIS WARRANTY SHALL NOT UNDER THE FOLLOWING CONDITIONS:**

- When the use, care and operation of the product has not been done in accordance with the instructions contained in the operating instructions.
- When the product has been used outside of its capacity, abused, beaten, exposed to moisture, ingressed by any liquid or corrosive substance, as well as any other fault attributable to the consumer.
- When the product has been disassembled, modified or repaired by unauthorized persons.
- The breakage and / or replacement of the bulb are not part of the warranty for this product.

In our continuous improvement process, Ortotecsa reserves the right to modify any characteristic without prior notice.